

CINCOM[®]

All For The Healthcare



Intelligence Eye Massager



Thank you for using the product. Please read the instructions carefully before using the product and follow the instructions to use it properly.

Catalogue

Thanks for purchasing our product, please read carefully before use. specially for attention part, and make sure of using correctly. Please keep it to check in further time.

Brief introduction	01
Graphene introduction	02
Infrared sensing	03
Application	04
Packing list	05
Use method	06
Function	07
Working mode & attentions	08
Basic fault & troubleshooting method	10
Maintenance	11
Ship&storage&specification	12
Contact us	13

If the design and changes of product functions are involved, without prior notice, please refer to the actual product.

Introduction

Product principle

The intelligent air pressure eye massage adopts the principle of traditional Chinese meridian and beauty, combined with the micro-chip control technology, improve eye blood circulation by soft warm, vibration and simulation of acupuncture on the eye. Built-in speaker and soothing music and relax you even more.

1. Intelligent air pressure : vibration to acupuncture point of eye and eyelid, Dredges the eye and head nervous system
2. Constant temperature soft warm: it takes by carbon fiber heating element, portable & constant temperature for eye massage, activates eye blood circulation, increases eye nutrient supply ability, effectively relieves eye dryness.
3. Vibration massage : vibrating the eye and eyelid, relieving blood circulation to reduce the elasticity of the ciliary muscle, enhancing the ability of the eye muscle to adjust.
4. Music: Relax with soothing music in 15 minutes.

Graphene

Graphene A planar film composed of carbon atoms. It is currently the world's lightest and most tough new nanomaterial due to its special nanostructure and high quality physicochemical properties. Graphene has important application prospects in electronics, optics, magnetism, biomedicine and drug delivery, and has been recognized as a "future revolution material"

Graphene Feature

portable convenient



stable performance





long life



Infrared portable physiotherapy:

Most close to body of wave physiotherapist:

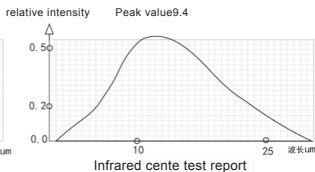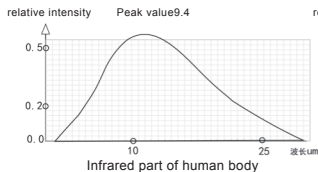


When graphene is heated, the far-infrared wavelength emitted is very close to the far-infrared wavelength of the human body, both of which are 6-14um. When the two wavelengths interact equally, resonance occurs. When resonance occurs with cells, The far-infrared heat energy penetrates into the subcutaneous tissue, uses capillary expansion, promotes blood circulation, strengthens tissue metabolism, increases regenerative capacity, and improves immunity of the body, thereby achieving medical care.

Light of life --Infrared sensing:

Light of life power of health:

Far-infrared rays have been discovered as early as 1800 years ago. It is an electromagnetic wave that is invisible to the naked eye. Part of the sun is a dizzying effect of life and growth. It is inextricably linked to the life of the Earth. Aerospace scientists have been in existence. The human living conditions in the spacecraft under vacuum, low temperature and overload conditions are investigated. It is known that far infrared rays with a wavelength of 6-12 microns in sunlight are essential for human survival. This wavelength of light is related to the human body. The far-infrared rays emitted are of the same wavelength, so that the human body can absorb the light energy and exchange and transmit the human body. Therefore, people call this long-range far infrared ray "the light of life."

Introduction of 10 types'far infrared physiotherapy:

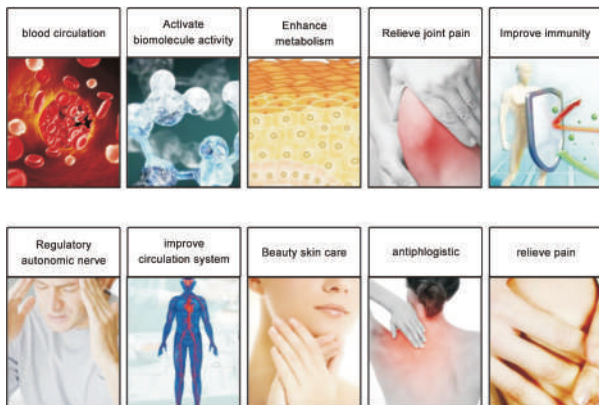

Application

- ▶ Large, middle and primary school students, myopia, hyperopia, over-eye or long-wearing glasses.
- ▶ Computer operators, accountants, designers, secretaries, teachers, long-time eye workers.
- ▶ Eye bags dark circles insomnia.



Large, middle and primary school students, myopia



hyperopia



over-eye or long-wearing glasses



Designer secretary Teacher



Driver



long-time eye workers



Computer operator



Insomnia



Long-term mobile phone player

Packing list



User manual
(including warranty/certificate)



USB Cable



Storage bag

Prepare before using

- ▶ Please confirm if the product has electricity. If the battery is too low and the machine function is turned off or cannot be turned on, please charge it before using it. Be sure to unplug the power adapter after charging, to avoid charging for a long time.

Charge

1. Connect the power adapter to the product USB port.
2. The charging LED light is red during charging, and the indicator light is blue when fully charged (or the LGD display battery power symbol is full).
3. During the charging process, the product does not support the booting work (charging in the working state, the product starts the self-protection program, and automatically stops working).
4. The main unit is powered by a built-in rechargeable battery, and the charging time is 3-4 hours.
5. Prompt for low battery charge during use. Please pause the use and recharge the battery.

Multiple charging methods

charging methods	computer USB port
	mobile
	standard power adaptor(recommend)


working time

The default working time for booting is 15 minutes. If you need to extend the usage time, you need to boot it again.

Function :

- 1.ON/OFF button: long press for 3 seconds to boot, long press for 3 seconds to shut down,Connect to Bluetooth,Double click to play music,Double-click ON/OFF again to pause playing music.Short press to switch modes.
2. The Bluetooth name is “Eye Massager”. Please connect Bluetooth before playing music for the first time. When Bluetooth is connected, speaker reminder: Connected; when Bluetooth is disconnected, speaker reminder: Disconnected.
3. 3D vibration mode-left side -right side- left&right side (upgrade motor)
4. This product cannot adjust the volume, so set the volume.80% of the volume, suitable for most people, other people Adjust the volume level with your phone.
5. With low battery voice reminder function,When the battery voltage is low, the speaker will prompt “Battery low”


6. Infrared operation:

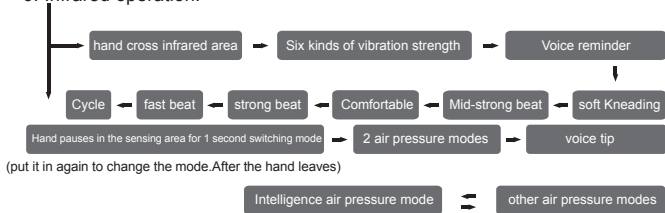

Attentions:

The six force function switches can only be switched in the vibration mode, otherwise it is invalid. Switching between the two air pressure modes only switches in the air pressure mode, otherwise it is invalid.

Function mode

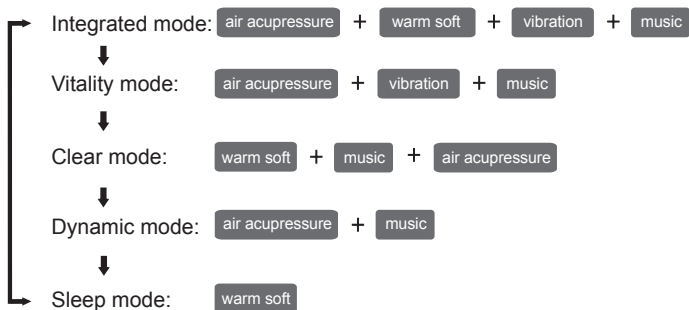

Attentions

To ensure that you use the massager correctly, please refer to the precautions before use:

1. Do not use the eye massager while walking outdoors or driving.
2. Do not use the eye massager in a hot or humid place.
3. There are symptoms of eye diseases that cannot be used. Major eye patients such as eye surgery, glaucoma, cataract, and retinal detachment have been banned.
4. Patients with deep stroke, hypertension and hypotension should be used with caution. Please consult a professional doctor. If there is any discomfort during use, please stop using it immediately.

5. Persons with mental and physical discomfort or mental retardation should be used with caution. Please use it under the supervision of the accompanying staff. If there is any discomfort during use, please stop using it immediately.
6. Please avoid using the eye massage device when drinking alcohol or physical discomfort.
7. Eye massager has a hot compress function, which must be taken care of when using temperature sensitive people.
8. Put the massage device on the eye area then turn it on again to avoid damage to the air bag.
9. The contact lenses, lenses, hairpins, hair accessories, etc. worn before use should be removed first.
10. This massage device is a personal health care product and cannot replace medicines and medical equipment.
11. If you feel pain or serious discomfort during use, please stop using it immediately and consult a physician.
12. This massager is designed as a home electronic health care product and is not intended for use in diagnostic, therapeutic and other commercial applications.
13. If cleaning is required, wipe it with a cotton cloth soaked with medical alcohol. Do not wash the product with water. Do not soak it in liquid.
14. Do not disassemble this product by yourself.
15. Please use the power adapter configured by our company and pay attention to electricity safety.
16. Please use this product for closed eyes. After use, the eyes will be blurred. After closing eyes for 1-3 minutes, the eyes will slowly clear.

Basic fault and troubleshooting method

basic fault	possible cause	troubleshooting method
failure to start up	too low power	charge in time
	Power failure	Please press the power button and hold for 5 seconds to turn on.
	unable to charge	damage adapter
low massage strength	too low power	charge in time
Sudden interruption of operation	The massager will shut off automatically in 15 minutes.	Press the power button to restart it if you need more massage.
Noise during use	<ol style="list-style-type: none"> 1.The sound of "buzzing" when the airbag is inflated 2. It is the normal sound of the machine when working. 3. It will not affect the use 	Don't need to do anything. If it is unusual, contact us via email: service@cincomhealth.com .

Maintenance

- 1.Clean with a dry soft cloth if there is dirt on the surface.
- 2.Clean with a soft cloth dampened with medical alcohol if particularly dirty.
- 3.Make the product in a safe place after using.
- 4.Take measures to prevent moisture, anti-rat,anti-pollution.
- 5.Please ask a professional technician to repair it if troubleshooting . Do not disassemble it at will.

Atmospheric pressure: 86 kpa~106 kpa

Ambient temperature: -10°C~65°C

Relative humidity: ≤80% RH, no corrosive gas and well ventilated room.

Transportation: This product is equipped with simple anti-shock facilities, suitable for air, rail, road and ship transportation. It should avoid rain, snow splash, inversion and collision. When the storage period is more than 6 months, it should be checked out of the box and energized for half an hour. It cannot be close to the ground, walls and roof. The room should be well ventilated to avoid strong sunlight and corrosive gas erosion.

Specification

Product Name: Eye Massager

Model NO.: CM-035E

Rating voltage: 3.7V

Rating power: 3.75W

Input voltage:5V/1A

Contact us

We provide 2-year warranty and lifelong after-sale service. If you have any problem during using of this product, please feel free to contact us. Our professional customer service team will help you!

E-mail: service@cincomhealth.com

Note: Please write the order number together with the problems you meet in the mail



 Cincom

 www.cincomhealth.com

Email: service@cincomhealth.com



All For The Healthcare

SHENZHEN CINCOM E-COMMERCE CO., LTD.

Address: B806, building 11, gongyuandadixiaoqu, Longxiangdadao
longchengjiedao, longgangqu, shenzhen, guangdong, China

E-mail: service@cincomhealth.com

Website: www.cincomhealth.com



RoHS
Made in China

